



Yes, you can.®

> Invacare® *Dragon Vertic* 

Access all areas – with the push on a button

The *Invacare® Dragon Vertic* is an agile power chair that not only provides mobility to an active user, but will bring you face to face with other people. Its comfortable seating system easily moves you into a standing position by a simple joystick command, providing increased access and excellent mobility.





➤ **Mobility when upright**

The *Invacare® Dragon Vertic* greatly enhances activity in your daily life with the vertical function. With mini front castors for extra stability you can even drive whilst in the upright position giving you increased access and mobility.

➤ **Verticalization**

The process of going upright is always from the seating position (driving) into vertical. The speed for moving up and down can be controlled via the joystick.

➤ **Easy transfer**

The knee fixation can be swivelled left or right, or be taken away entirely for easier transfer.

The transfer is also facilitated by swivelling-up arm rests and foot plates.

➤ **Comfortable positioning**

Seat depth, leg rest length, back rest tension and height can be adapted to the needs of the user, and the angle of the back rest can be electrically adjusted for optimum comfort when driving and when upright.



Features and options



The arm rest can be turned upward for an easier transfer.



The head rest is individually adjustable in width, height and horizontal position.



The sliding back rest minimises the shearing forces. Its height is adjustable.



The hip support is width adjustable and removable for transfer.



The adjustable back rest angle optimises comfort in seating and upright position.



The foot plates can be swivelled up for easy transfer.



Mini castors allow driving in the upright position.



Puncture proof rear tyres available.

Features and options



Tray

Convenient tray for work and for trunk support, can be swivelled down.



Knee support

The knee support is adjustable in depth, width and height. It can be swivelled left and right, or can be taken off entirely.



Shark II


Easy steering with **Shark II**. Variable speed for vertical function and for back rest angle adjustment.



On the road version

Front and rear lights provide safety when driving on public roads. The light is activated via the Shark controller.

Technical data

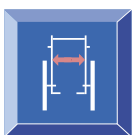

38,5 - 46 cm


45 cm incl. armpads



38 - 46 cm



50 / 55 cm

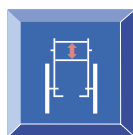


60 / 65 cm



26 - 35 cm




31 - 45 cm

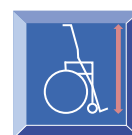

167 / 172 cm



6° - 22°

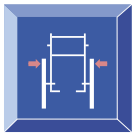

48,5 - 56 cm*1

62,5 cm incl. armpads




116 / 121 cm



105 cm*2



85 cm

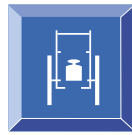


131 kg



100 kg



2 x 12V 60 Ah

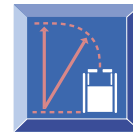

83 cm



4 cm




17,6 %



2 x 180 W



25 km*3

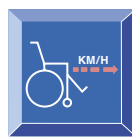

6 km/h

*1 61 cm for driving base
*2 incl. leg rest and antitippers
*3 acc. to ISO 7176

Frame colour



Blue metallic

All measurement have been done under the Invacare standard test conditions, unless quoted differently.
All information is believed to be correct at the time of print. Invacare® reserves the right to alter product specifications without prior consultation.



INVACARE® Deutschland GmbH
Kleiststraße 49 · 32457 Porta Westfalica
Telefon + 49 (0) 57 31 / 7 54-0 · Telefax + 49 (0) 57 31 / 7 54-5 21 91
E-Mail info@invacare.de · www.invacare.de

