



*Argon*

First time, second time, everytime



QUICKIE

INNOVATE.  
BUILD FOR THE FUTURE

# Form follows function

A new era has begun: Argon, the

The Argon from Quickie treads new paths with regards to flexibility and adjustability. The Argon concept is based on the maximum in individuality with the minimum time and expense needed to make the adjustment. In addition, the Argon will impress you with its extremely light weight, excellent drive characteristics and a powerful look. The new Argon – ideal for all the requirements of daily life, during leisure time, at work and playing sport.



# QUICKIE

forward-looking wheelchair concept.

## Take your seat! The ultimate in adjustability.

Good drive characteristics almost go without saying on a fixed frame wheelchair. If you can also carry out individual adjustments and positioning in next to no time, then no desire is left unfulfilled.

The individual adjustment of the back can be carried out while you are sitting in the wheelchair – in this way, you can feel every change immediately. So the user can find his own optimum and most comfortable seating position without having to keep getting in and out, in the shortest possible time.



1



2



3



4

- 1-3 Opening the backrest clamp
- 4 Adjustment of backrest angle

The centre of gravity of the new Argon is infinitely adjustable. The back height is easy to adjust, and the seat depth or front seat height can also be adjusted to meet individual requirements. The seat angle can be altered in 2° steps, simply by turning the castor fork.



1



2



3



4

- 1 Centre of gravity
- 2 Fixed castor housing
- 3 Height-adjustable castor housing
- 4 Rear seat height adjustment

As an option, the Argon is available in various heights, widths, with fixed or adjustable castors. A multitude of option possibilities ensure that the new Argon is the perfect chair for every situation – a true all-rounder.

### You show style

The design of the new Argon sets new standards. The curved lines give this active chair a light and open overall look. With 27 different frame colours through to carbon look sideguards all combine to offer just the right solution for every taste. Even the newly-designed castor wheel receiver emphasizes the universal design of the new Argon – Form and Function in perfect harmony.




# Argon

## The therapist and Argon – a perfect team

The optimum adjustment of a wheelchair to meet the requirements of the user often causes the therapist problems. Compromises because of the chair construction or limited adjustment possibilities are, sadly, the order of the day. Perfect seat positioning or perfect transfer of power can usually only be achieved through costly discussions and time-intensive adjustment sittings.



This is completely different with the new Argon. The Argon sets standards in relation to easy adjustment, flexibility and thereby the optimum for individuality. From initial supply to individual adjustment of an active chair, the new Argon can be easily and quickly adjusted to the requirements of the user.

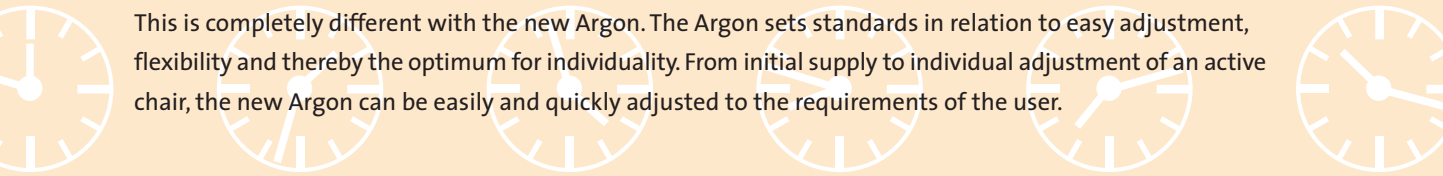

Through the completely new overall concept of the chair, the Argon can be adjusted to the user, without having to constantly get into and out of the chair. In this way, the therapist gets constant feedback and can find the ideal position for the user with the minimum amount of alteration – conceivably simple and rapid, too.



### **Pure performance**

The fixed caster connection ensures minimum weight with optimum function. With the height-adjustable caster receiver which is available as an option, the Argon meets all therapeutic requirements.

### **The frame – as individual as the concept**

Also when it's a question of frames, the Argon is geared completely to individuality. Whether its with short or long frame angle (95 and 105) or the V-shaped frame curved inwards (0, 2 and 4 cm), the main priority is always the optimum adjustment to the meet the requirements of the user.

# QUICKIE



## Overview of the new flexibility:

|                                    |                                      |
|------------------------------------|--------------------------------------|
| Seat width                         | from 30 – 50 cm                      |
| Seat depth                         | from 30 – 50 cm                      |
| Back height                        | from 30 – 45 cm +/- 2,5 cm           |
| 27 individual frame colours        |                                      |
| 6 individual colours of upholstery |                                      |
| 2 frame angles                     | 95°, 105°                            |
| 3 frame inset                      | 0 cm, 2 cm, 4 cm                     |
| 2 versions of frame                | fixed and variable front seat height |
| Weight                             | from 9,3 kg on                       |
| Optional seat depth growth         | 2,5 and 5 cm                         |
| Suspension system                  |                                      |
| 5 different footrest               |                                      |

New lifestyle options available:

26" Lightweight wheel, 24" Offroad wheels, 4" and 5" design front wheels

**STPS**  
Sunrise Training and Education Programmes

**KNOWLEDGE FOR THE FUTURE**

**SUNRISE  
MEDICAL**

Sunrise Medical Limited  
High Street Wollaston  
West Midlands DY8 4PS  
England

Phone +44 (0) 1384 446622

Fax +44 (0) 1384 446644

[www.sunrisemedical.com](http://www.sunrisemedical.com)



**INNOVATE.  
BUILD FOR THE FUTURE**