




Proudly made
and manufactured
in Canada!

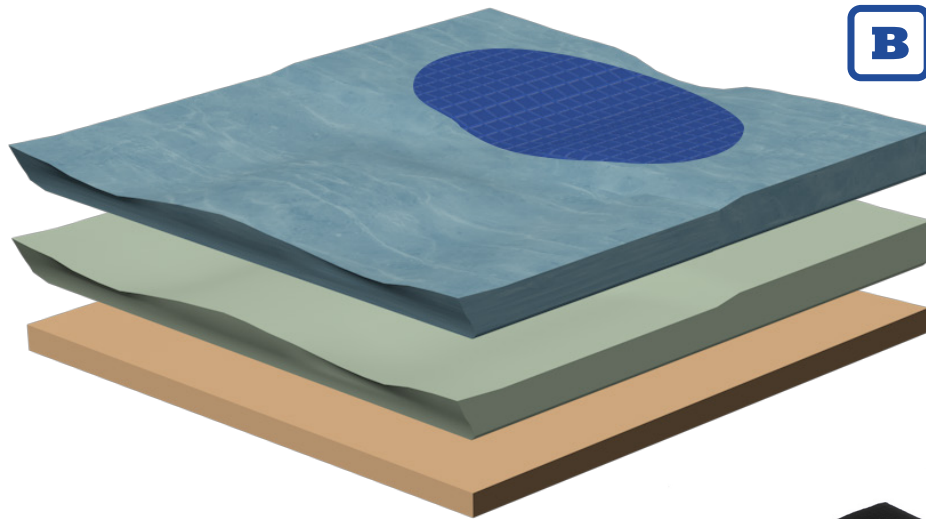

Product Specifications

CONSTRUCTION

Closed cell foam base and lateral stabilizing wedges. High resiliency foam for the base, and ischial shelf. Core is comprised of a contoured injected foam. The top layer features a gel infused memory foam at the front of the cushion and a zoned Pillow Top Geo-Matrix™ gel insert along the back half of the cushion.

COVERS

Two outer cover options: Breathable Spacetex fabric with incontinent protection OR Elastatex soft touch outer incontinent cover.

Inner cover: Elastatex breathable, wipeable, antibacterial, fluid proof cover.

STANDARD SIZES

35x35, 35x38, 35x40, 35x43, 35x48, 35x50, 38x35, 38x38, 38x40, 38x43, 38x45, 38x48, 38x50, 40x35, 40x38, 40x40, 40x43, 40x45, 40x48, 40x50, 40x53, 40x55, 43x35, 43x38, 43x40, 43x43, 43x45, 43x48, 43x50, 43x55, 45x30, 45x33, 54x35, 45x38, 45x40, 45x43, 45x45, 45x48, 45x50, 45x53, 45x55, 48x35, 48x38, 48x40, 48x43, 48x45, 48x48, 48x50, 48x53, 48x55, 50x35, 50x38, 50x40, 50x43, 50x45, 50x48, 50x50, 50x53, 50x55, 53x35, 53x38, 53x40, 53x43, 53x45, 53x48, 53x50, 53x53, 53x55, 55x35, 55x38, 55x40, 55x43, 55x45, 55x48, 55x50, 55x53, 55x55, 58x40, 58x43, 58x45, 58x48, 58x50, 58x53, 58x55, 60x40, 60x43, 60x45, 60x48, 60x50, 60x53, 60x55, 60x55

CUSTOM SIZES

All dimensions available

PRODUCT CODE

GV(+size) eg. GV4545

WEIGHT CAPACITY

Recommended Maximum 125kg


The Geo-Matrix™ Vinga cushion perfectly combines the pressure and shear reduction of Geo-Matrix™ gel with the plush comfort of gel infused memory foam.

The Vinga was designed to maximize comfort and protection while maintaining a stable posture. With its multi foam core and the firm ischial shelf the Vinga finds the balance between comfort and stability.

The Vinga is designed to be light-weight and maximize comfort, while maintaining stable posture.

Each Vinga cushion comes complete with either a breathable or soft touch incontinent cover, as well as a 4 way stretch protective inner cover.

Blake Medical esteemed distributor:

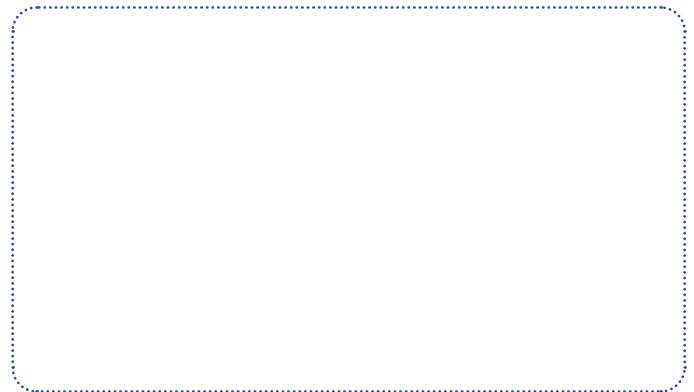

FOAM CUSHION
Geo-Matrix Vinga




WASHING INSTRUCTIONS

Covers: 60°C in washing machine. Hang to air dry. Do not bleach.

Geo-Matrix™ gel and foam: Wipe down gel using any cleaner that does not contain bleach and allow to air dry.

For more specific cleaning instructions please visit our website or user's manual.

All of Blake Medical products are free from phthalates and latex. This includes all medical fabrics, Geo-Matrix™ gel, moulded polyurethane, standard polyurethane and specialty foam. All cushions meet EN-1021 part 1 & 2 requirements for fire retardancy.

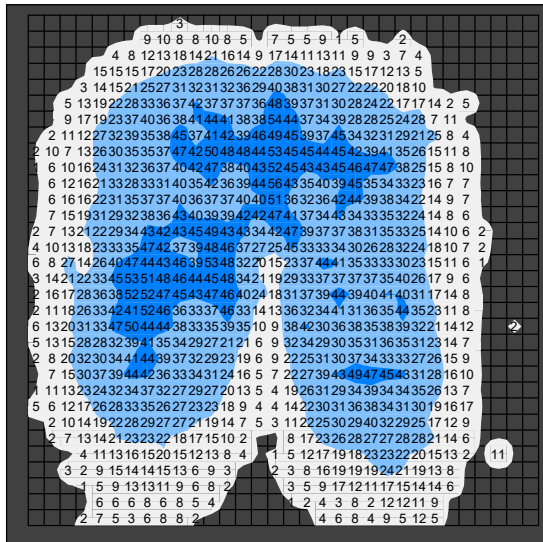

Client: Female 53kg

ADDITIONAL BENEFITS

- Cushion remains 2-3 degrees cooler under seated body temperature to assist with microclimate management.
- Minimizes shear through column segmentation.
- Geo-Matrix™ Gel does not change structure in climates of extreme heat nor extreme cold.
- Geo-Matrix™ gel is recyclable through Blake Medical.

FOAM CUSHION
Geo-Matrix Vinga

Advantages of Geo-Matrix™ Gel

DURABLE	COMFORTABLE	DISINFECT-ABLE	FRICTION & SHEAR REDUCING	INCREASES SIT TIMES
MICROCLIMATE MANAGEMENT	PRESSURE REDISTRIBUTING	VIBRATION DAMPENING	WON'T PUNCTURE	

