

Drive
Medical
Specialist
Seating





Seating & Backs



SHELL SIZE	WEIGHT LIMITS
30.5cm - 50.8cm (12" - 20")	113kg (17.5st)
55.9cm - 61cm (22" - 24")	295kg (46st) (with reinforcement bars)



WIDTH	DEPTH			
	30.5cm (12")	35.6cm (14")	40.6cm (16")	45.7cm (18")
30.5cm (12")	PFDSV12I2BE	PFDSV12I4BE		
35.6cm (14")		PFDSV14I4BE	PFDSV14I6BE	
40.6cm (16")		PFDSV16I4BE	PFDSV16I6BE	
45.7cm (18")			PFDSV18I6BE	PFDSV18I8BE
50.8cm (20")			PFDSV20I6BE	PFDSV20I8BE

Remove (V) from the product code should you not require VERSALocks with the Drop Seat



1. VERSATEC BACK SYSTEM

- » Shell is fabricated from highly durable, premium grade Acrylic/PVC
- » Partially pre-drilled holes maintain a clean look from the outside and can then be easily punched out for precise placement of mounting hardware component parts
- » Available in either Standard or Deep Contour for enhanced trunk stability
- » Deep contour panels can be mounted on either the Standard or Tall Shell Sizes
- » Installs in seconds with Tension Hook hardware
- » System has been fully crash tested in compliance with ISO16840 Part 4 DIS, incorporating the uniquely patented VERSALock system
- » Utilises hook and bracket combinations to ensure a greater range of forward and drop-angle adjustment to accommodate growth and individual needs
- » Side to side adjustments for asymmetrical positioning
- » Fewer back sizes accommodate multiple chair widths
- » Optional VERSALocks provide an easy and certain way to mark your back location when removed from the frame

2. BACK CUSHION

- » The contoured back cushion is covered in a two-way stretch, vapour permeable, water resistant fabric manufactured in the UK
- » The network of channels allows air to circulate helping to reduce the build up of heat and moisture
- » Velcro® strips make it easy to remove for fast and easy cleaning
- » Lumbar and lateral inserts included
- » For sizes please see order guide & size information at back of the brochure

3. ADJUSTABLE DROP SEAT

- » Easy to use, solid, fully adjustable seat base provides a firm base of support for wheelchair cushions
- » Adjustable seat width accommodates growth and stability
- » Angle adjustable through Tension Hook hardware
- » Supplied with VERSALocks as standard

4. TENSION HOOK HARDWARE

- » Hook and loop to aid cushion positioning
- » Available in 10 sizes

Headrests & Brackets

1. SOFT HEADREST

- » Designed for both comfort and support
- » Removable and washable cover
- » Low height ensures discreetness for user
- » Cost effective pad and bracket combinations

1



HEIGHT	WIDTH	HEADREST W/STD. REMOVABLE BRACKET	HEADREST W/ POSALINC REMOVABLE BRACKET
12.7cm (5")	28cm (11")	P1802162BE	P1802262BE
14cm (5.5")	35.5cm (14")	P1802462BE	P1802562BE

2. TRI-PAD HEADREST

- » Provides more aggressive head support without affecting sensory perception
- » Curved pad with occipital contour at the bottom of the middle pad
- » Optional forehead strap

2



CODE	SIZE	HEIGHT	WIDTH	DEPTH
P1808762BE	Regular	Side Pads: 3.8cm-7.6cm (1.5"-3") Middle Pad: 10.2cm (4")	Middle Pad: 11.4cm (4.5")	Side Pads: 12.7cm (5")
P1809262BE	Large	Side Pads: 4.5cm-8.9cm (1.75"-3.5") Middle Pad: 15.2cm (6")	Middle Pad: 14cm (5.5")	Side Pads: 15.2cm (6")

3. POSALINC ADJUSTABLE HEADREST BRACKET

- » Full range of adjustment to position the headrest pad in precise location
- » Removable brackets
- » Can accommodate a range of headrest pads
- » 45° Vertical Rod

3



CODE	PRODUCT
P18191BE	POSALinc Bracket

4



CODE	BAR	POST
P18053BE	Short	Long
P18051BE	Long	Long

Laterals

1



CODE	HEIGHT	SLOT SPACING	LENGTH
P15224BE	5.1cm (2")	2.5cm (1")	10.1cm (4")
P15225BE			12.7cm (5")
P15226BE			15.2cm (6")
P15228BE			20.3cm (8")

2




CODE	SIZE
P2516062LBE (Left) P2516062RBE (Right)	To fit shell size 30.5cm - 35.5cm (12" - 14")
P2516262LBE (Left) P2516262RBE (Right)	To fit shell size 40.5cm - 51cm (16" - 20")

3



CODE	SIZE
P2525062LBE (Left) P2525062RBE (Right)	To fit shell size 30.5cm - 35.5cm (12" - 14")
P2525262LBE (Left) P2525262RBE (Right)	To fit shell size 40.5cm - 45.7cm (16" - 18")

4

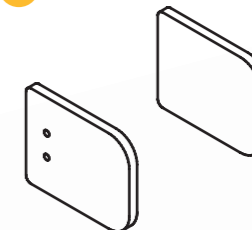

CODE	HEIGHT	DEPTH
23257BE	10.1cm (4")	15.2cm (6")
P2640862BE	10.1cm (4")	20.3cm (8")

5



CODE	HEIGHT	SLOT SPACING	MOUNTING
P15114LBE (Left) P15114RBE (Right)	12.7cm (5")	2.5cm (1")	Straight
P15117LBE (Left) P15117RBE (Right)	12.7cm (5")	2.5cm (1")	1" Offset

1. STANDARD LATERAL BRACKETS

- » Designed to mount pads off solid backs or seats
- » Available in a range of lengths to accommodate varieties of back and seat thicknesses
- » Slots allow for both side to side and forward to back adjustment
- » 90° angle for easy mounting to solid flat backs

2. FIXED LATERAL BRACKETS

- » Designed to mount lateral pads off wheelchair backs
- » Available as pairs or single units, in two sizes

3. QUICKPLUS SWING-AWAY LATERAL BRACKETS

- » Help facilitate quick and easy transfers with side to side lateral adjustment
- » Steel manufactured brackets with black powder coated finish
- » Angle adjustment to rotate front of pad in or out
- » Available as pairs or single units, in two sizes

4. FLAT HIP PADS

- » Strong yet flexible offering increased comfort, support, and adjustability


5. REMOVABLE HIP PAD BRACKETS

- » Slim style design with easy one-hand release spring-loaded locking mechanism

Belts & Pelvic Positioners

1. DYNAFORM POSTURE SUPPORT



- » Comfortable support whilst offering an adjustable level of dynamic positioning
- » Front panel is removable for laundering
- » Includes side release buckles and belt anchors for easy adjustment and attachment to a solid back
- » Lycra® and nylon covered neoprene conforms to the shape of the body whilst allowing freedom of movement

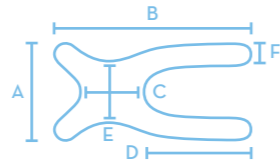

CODE	SIZE	A	B	C	D	E	F
P14850BE	XX Small	12.7cm (5")	25.4cm (10")	10.2cm (4")	12.7cm (5")	3.2cm (1.25")	4.4cm (1.75")
P14851BE	Extra Small	20.3cm (8")	35.6cm (14")	10.2cm (4")	21.6cm (8.5")	10.2cm (4")	5.1cm (2")
P14852BE	Small	22.8cm (9")	43.2cm (17")	12.7cm (5")	25.4cm (10")	12.7cm (5")	5.1cm (2")
P14853BE	Medium	25.4cm (10")	50.8cm (20")	15.2cm (6")	29.2cm (11.5")	15.2cm (6")	5.1cm (2")
P14854BE	Large	27.9cm (11")	55.9cm (22")	17.7cm (7")	30.5cm (12")	17.7cm (7")	6.3cm (2.5")

2. SLIM-CUT DYNAFORM POSTURE SUPPORT


- » Designed for female clients or those who want less material covering the chest area
- » Lycra® and nylon covered neoprene conforms to the shape of the body whilst allowing freedom of movement

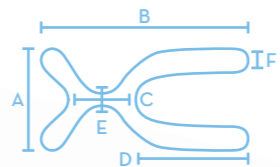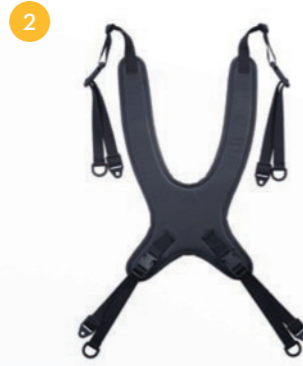


CODE	SIZE	A	B	C	D	E	F
P14856BE	Extra Small	20.3cm (8")	35.6cm (14")	10.2cm (4")	21.6cm (8.5")	2.5cm (1")	5.1cm (2")
P14857BE	Small	22.8cm (9")	43.2cm (17")	12.7cm (5")	25.4cm (10")	2.5cm (1")	5.1cm (2")
P14858BE	Medium	25.4cm (10")	50.8cm (20")	15.2cm (6")	29.2cm (11.5")	3.2cm (1.25")	5.1cm (2")
P14859BE	Large	27.9cm (11")	55.9cm (22")	17.7cm (7")	30.5cm (12")	3.2cm (1.25")	6.3cm (2.5")

3. UPPER BODY DYNAFORM POSTURE SUPPORT

- » Designed to eliminate obstruction in the lower trunk area for feeder tubes, hip belts, support pads etc
- » Centre buckle design offers quick release for transferring the client as well as easy width adjustment
- » Lycra® and nylon covered neoprene front panel
- » Includes centre and side release buckles with belt anchors for easy adjustment and attachment to a solid back

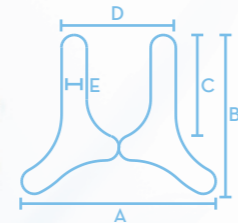


CODE	SIZE	A	B	C	D	E
P14867BE	Small	35.6cm (14.25")	26.7cm (10.5")	15.9cm (6.25")	21.6cm (8.5")	5.1cm (2")
P14868BE	Medium	45.7cm (18")	34.9cm (13.25")	20.3cm (8")	26cm (10.25")	5.7cm (2.25")
P14869BE	Large	48.9cm (19.25")	40cm (15.75")	25.4cm (10")	28.6cm (11.25")	5.7cm (2.25")

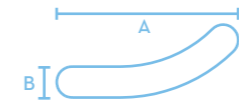


CODE	HARNESS SIZE	WEBBING SIZE	A	B
P14872BE	Small	2.5cm (1")	27.9cm (11")	5.1cm (2")
P14873BE	Medium	2.5cm (1")	34.3cm (13.5")	5.7cm (2.25")
P14874BE	Large	2.5cm (1")	40.6cm (16")	5.7cm (2.25")



CODE	SIZE	HEIGHT	WIDTH	RECOMMENDED CHEST WIDTH
P14927BE	Medium	10.2cm (4")	50.8cm (20")	30.5 - 35.6cm (12 - 14")
P14928BE	Large	10.2cm (4")	55.9cm (22")	35.6 - 40.6cm (14 - 16")

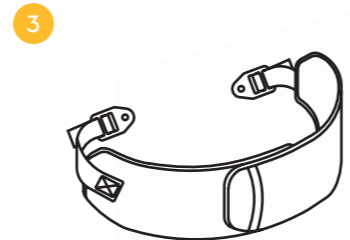

CODE	SIZE	HEIGHT	RECOMMENDED CHEST WIDTH
P14942BE	Medium	10.2cm (4")	25.4 - 40.6cm (10 - 16")
P14943BE	Large	12.7cm (5")	30.5 - 45.7cm (12 - 18")

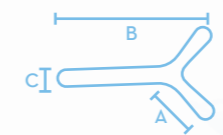

CODE	SIZE	A	B	C
P14862BE	Small	17.7cm (7")	45.7cm (18")	5.1cm (2")
P14863BE	Medium	20.3cm (8")	50.8cm (20")	5.1cm (2")
P14864BE	Large	25.4cm (10")	60.9cm (24")	6.3cm (2.5")

1. BACKPACK STYLE DYNAFORM POSTURE SUPPORT

- » Designed to pull the shoulder girdle back and position the upper extremities from internal rotation
- » Offers alternative to standard style posture supports that can interfere with breast tissue
- » Lycra® and nylon covered neoprene front panel
- » Includes centre and side release buckles with belt anchors for easy adjustment and attachment to a solid back

2. DYNAFORM CHEST SUPPORT: ONE PIECE BUCKLE

- » Designed to prevent the client's trunk from falling forward, while stretching to allow movement and easier breathing
- » Maintains contact with the front and sides of the client and mounts to a solid back or to the wheelchair uprights
- » Provides the client with rotation of the shoulders and greater ability to reach forward
- » Easily adjustable and removable
- » Lycra® and nylon covered neoprene material stretches for dynamic movement

3. DYNAFORM CHEST SUPPORT: TWO PIECE HOOK/LOOP

- » Maintains contact with the front and sides of the client and mounts to a solid back or to the wheelchair uprights
- » Easily adjustable and removable
- » Lycra® and nylon covered neoprene material stretches for dynamic movement

4. Y-STYLE POSTURE SUPPORT

- » Designed to support and position for spinal/shoulder rotation requiring the pull and control on one body side only
- » Allows for close, comfortable fit around the shoulder and across the body
- » Apply to left or right side of the client, suitable for both male and female body contours

Belts & Pelvic Positioners

1. ARCUfit

The ARCUfit Pelvic Positioner incorporates a semi-rigid design that combines the positive features of lateral hip supports, padded hip belts, and sub-ASIS bars into one great product to position the pelvis.

- » Semi-rigid/flexible premium grade composite base combined with spring steel
- » Comfortable contoured foam
- » Adjustable width and buckle placement
- » Angle adjustable positioning
- » Choice of attaching brackets mount off solid seat/back or to the frame
- » Controls pelvic movement similar to a sub-ASIS bar but with greater comfort and provides lateral hip support
- » Provides consistent placement of support to maintain proper pelvic positioning
- » Spring action materials hold the positioner and buckles out of the way of wheels during transfers
- » Client and carers can easily access the ends of the positioners to fasten together

2. OFFSET-IN LATERAL BRACKETS


- » Designed to mount pads off solid backs or seats for clients who are narrower than the back or seat

1



CODE	SIZE	SEATED MEASUREMENT	BUCKLE TYPE
PI4951BE	Extra Small	38.1 - 48.3cm (15 - 19")	Push Button
PI4952BE	Small	48.3 - 60.9cm (19 - 24")	Push Button
PI4953BE	Medium	58.4 - 71.1cm (23 - 28")	Push Button
PI4954BE	Large	58.6 - 83.8cm (27 - 33")	Push Button

2

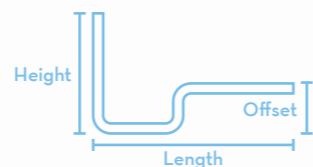

CODE	HEIGHT	SLOT SPACING	OFFSET	LENGTH
PI5820BE	5.1cm (2")	2.5cm (1")	1.25cm (0.5")	10.1cm (4")
PI5821BE	5.1cm (2")	2.5cm (1")	1.25cm (0.5")	12.7cm (5")
PI5822BE	5.1cm (2")	2.5cm (1")	2.5cm (1")	10.1cm (4")
PI5823BE	5.1cm (2")	2.5cm (1")	2.5cm (1")	12.7cm (5")

Please see page 5 for Lateral Bracket options

1



SINGLE PULL CODE	BUCKLE TYPE	SIZE	WEBBING SIZE	PAD LENGTH
PI4603BE	Side Release	Extra Small	2.5cm (1")	15.2cm (6")
PI4613BE		Small	3.8cm (1.5")	20.3cm (8")
PI4614BE		Medium	3.8cm (1.5")	25.4cm (10")
PI4624BE		Large	5.1cm (2")	30.5cm (12")
PI4633BE	Push Button	Extra Small	2.5cm (1")	15.2cm (6")
PI4643BE		Small	3.8cm (1.5") Spring	20.3cm (8")
PI4644BE		Medium	3.8cm (1.5") Spring	25.4cm (10")
PI4663BE		Small	3.8cm (1.5") Non-Spring	20.3cm (8")
PI4654BE		Large	5.1cm (2")	30.5cm (12")
PI4683BE	Lever Release	Large	5.1cm (2")	25.4cm (10")

2



CODE	BUCKLE TYPE	BELT SIZE	WEBBING SIZE	PAD LENGTH
PI4749BE	Push Button	Small	2.5cm (1")	13.9cm (5.5")
PI4750BE		Medium	3.8cm (1.5")	19.1cm (7.5")
PI4751BE		Large	5.1cm (2")	23.5cm (9.25")
PI4754BE	Plastic Side Release	Small	2.5cm (1")	13.9cm (5.5")
PI4755BE		Medium	3.8cm (1.5")	19.1cm (7.5")
PI4756BE		Large	5.1cm (2")	23.5cm (9.25")

1. HIP BELTS

Used to hold the pelvis in the correct position within the seating system, the padded hip belts have fabric and foam laminated pads cut and trimmed with a tapered design for exceptional comfort and streamlined appearance.

- » Choose from Single Pull or Dual Pull (pull forward)
- » Includes anchors attached for easy adjustment and attachment to the frame

2. HIP STABILISING BELT

- » Maintains contact with the front and sides of the client and mounts to a solid back or to the wheelchair
- » Easily adjustable and removable
- » Lycra® and nylon covered neoprene material stretches for dynamic movement

ABDUCTORS & ADDUCTORS


1. ABDUCTOR PADS

- » Designed to position the legs out of adduction and into better alignment
- » Can help improve weight bearing balance and stability
- » Available in either oval or wedge shape

2. PUSH BUTTON ABDUCTOR BRACKET

- » Designed for function and style, the flip away bracket flips the abductor underneath the front edge of the seat to allow for easy transfers
- » Low profile

1

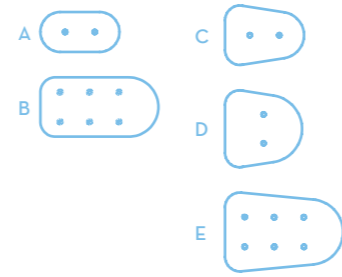

OVAL PAD					
CODE	SHAPE	WIDTH	LENGTH	HEIGHT	SUGGESTED TOP PLATE
P2224362BE	A	5.1cm (2")	10.2cm (4")	7.6cm (3")	Child
P2336362BE	B	7.6cm (3")	15.2cm (6")	7.6cm (3")	Long
WEDGE PAD					
CODE	SHAPE	WIDTH	LENGTH	HEIGHT	SUGGESTED TOP PLATE
P2234362BE	C	5.1 - 7.6cm (2 - 3")	10.2cm (4")	7.6cm (3")	Child
P2344362BE	D	7.6 - 10.2cm (3 - 4")	10.2cm (4")	7.6cm (3")	Short
P2346362BE	E	7.6 - 10.2cm (3 - 4")	15.2cm (6")	7.6cm (3")	Long

2



Child Top Plate Short Top Plate Long Top Plate

SIZE	SEAT THICKNESS	CHILD TOP PLATE	SHORT TOP PLATE	LONG TOP PLATE
Large	8.9 - 11.4cm (3.5 - 4.5")	P14343BE	P14313BE	P14323BE

1



CODE	TUBE DIAMETER
P13505LBE (Left) P13505RBE (Right)	2.2cm (0.875")
P13515LBE (Left) P13515RBE (Right)	2.5cm (1")

2



CODE	TUBE DIAMETER
P13405LBE (Left) P13405RBE (Right)	2.2cm (0.875")
P13415LBE (Left) P13415RBE (Right)	2.5cm (1")

1. POSAlinc KNEE ADDUCTORS

- » Adjustable design provides a pad and bracket combination that can be positioned out at the front of the knee for optimum support
- » Infinite adjustments allows optimum positioning
- » Can be mounted onto either the straight or curved tubing of the leg rest

2. BASIC KNEE ADDUCTORS

- » Cost effective option
- » Multi-adjustment pad positioning

ARM & FOOT ACCESSORIES

1. CONTOURED ARM SUPPORT

- » Durable lightweight base can be heat modified for additional contour
- » Comfortable nylon covered neoprene padding attached with Velcro® removes easily for washing
- » Pre-fabricated slots for wrist-strap option (wrist strap sold separately)
- » Swivel, tilt and removable mounting option
- » Includes fixed mounting plates to attach arm support securely onto wheelchair arm tube
- » Sized from joystick use to full length

2. POSAlinc MOUNTING HARDWARE

- » Swivel and tilt mounting
- » Ability to remove arm support with simple twist of ratchet style lever

1



CODE	LENGTH
PI9020BE	25.4cm (10")
PI9022BE	30.5cm (12")
PI9024BE	35.6cm (14")
PI9026BE	40.6cm (16")
PI9028BE	45.7cm (18")

2



CODE	LENGTH
PI9047BE	2.2cm (0.875")
PI9048BE	2.5cm (1")

Optional Padded Wrist Strap: PI9043

1



CODE	SIZE	FOOT BASE DIMENSION	WIDTH BETWEEN SLOTS
PI3161BE	Small	11.4 - 19.1cm (4.5 - 7.5")	6.3 - 8.9cm (2.5 - 3.5")
PI3162BE	Medium	12.7 - 22.9cm (5 - 9")	7.6 - 10.2cm (3 - 4")
PI3163BE	Large	13.9 - 25.4cm (5.5 - 10")	8.9 - 11.4cm (3.5 - 4.5")

2



CODE	SIZE	TOE	HEEL	LENGTH
22119ABE	Small	6.9cm (2.75")	6.9cm (2.75")	15.2cm (6")
PI3012BE	Medium	10.2cm (4")	7.6cm (3")	20.3cm (8")
PI3013BE	Large	11.4cm (4.5")	8.9cm (3.5")	25.4cm (10")

3



CODE	PANEL SIZE	WEBBING SIZE	PADWIDTH	PAD HEIGHT	FRAME SIZE
PI3175BE	Small	2.5cm (1")	48.3cm (19")	12.7cm (5")	30.5 - 38.1cm (12 - 15")
PI3176BE	Medium	2.5cm (1")	58.4cm (23")	16.5cm (6.5")	40.6 - 45.7cm (16 - 18")
PI3177BE	Large	3.8cm (1.5")	66.1cm (26")	20.3cm (8")	48.3 - 55.9cm (19 - 22")

1. DYNAFORM FOOT & ANKLE POSITIONERS

- » Allow dynamic stretch and movement, easily adjusting to a variety of footwear
- » Unique 3 point ankle attachment forms a cup for the heel whilst the side straps wrap around the ankle
- » Rigid plastic base is designed with multiple slots for adjusting the width
- » Contoured side release buckle for simple release
- » Ankle Positioner allows the foot to flex and rotate

2. SHOE HOLDERS

- » Holds the foot securely in position with Velcro® strap adjustment
- » Pads included for added comfort whilst easy to remove for cleaning
- » Contoured side release buckles conform to the shape of the foot and provide a simple release

3. CALF PANELS

- » Designed to keep legs in place
- » Comfortable Neoprene panel with easy to adjust Velcro® straps to fit most chair widths
- » Available in 3 sizes

HARDWARE

1. VERSALock SECURING CLIPS






- » Unique, patented design provides an easy to use yet secure method to lock a solid seat or back onto a wheelchair frame
- » The top portion of the VERSALock slides over the drop hook hardware with a spring-loaded snap-action to ensure it is locked into position
- » Made of durable injection moulded plastic
- » Easy snap-on attachment

1



CODE	QUANTITY
P8001002BE	2 pack
P8001004BE	4 pack

ORDERING GUIDE

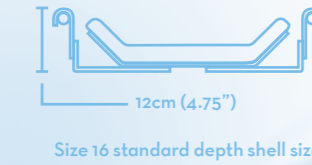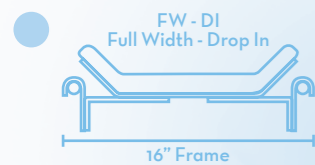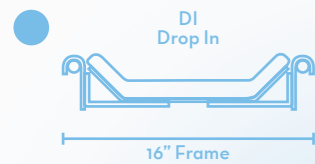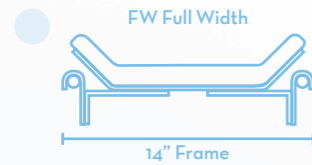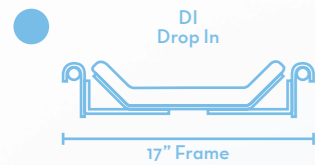


Fewer back sizes accommodate multiple chair widths

SCALE-SHELL SIZE	30.5CM (12")	35.6CM (14")	40.6CM (16")	45.7CM (18")	50.8CM (20")	55.8CM (22")	61CM (24")
30.5cm (12")	DI	FW					
33cm (13")	DI	FW - DI					
35.6cm (14")		DI	FW				
38.1cm (15")		DI	FW - DI				
40.6cm (16")			DI	FW			
43.2cm (17")			DI	FW - DI			
45.7cm (18")				DI	FW		
48.3cm (19")				DI	FW - DI		
50.8cm (20")					DI	FW	
53.3cm (21")					DI	FW - DI	
55.9cm (22")						DI	FW
58.4cm (23")						DI	FW - DI
61cm (24")							DI

SIZE INFORMATION

FRAME RANGE	SHELL SIZE	STANDARD	TALL	DEEP CONTOUR	
		STANDARD BACK W/ VERSALOCK	TALL BACK W/ VERSALOCK	STANDARD DEEP CONTOUR W/ VERSALOCK	TALL DEEP CONTOUR W/ VERSALOCK
30.5 - 33cm (12 - 13")	30.5cm (12")	PSV12BE	PTV12BE	PSV12DCBE	PTV12DCBE
30.5 - 38.1cm (12 - 15")	35.6cm (14")	PSV14BE	PTV14BE	PSV14DCBE	PTV14DCBE
35.6 - 43.2cm (14 - 17")	40.6cm (16")	PSV16BE	PTV16BE	PSV16DCBE	PTV16DCBE
40.4 - 48.3cm (16 - 19")	45.7cm (18")	PSV18BE	PTV18BE	PSV18DCBE	PTV18DCBE
45.7 - 53.3cm (18 - 21")	50.8cm (20")	PSV20BE	PTV20BE	PSV20DCBE	PTV20DCBE
50.8 - 58.4cm (20 - 23")	55.8cm (22")	PSV22BE		PSV22DCBE	
55.9 - 61cm (22 - 24")	61cm (24")	PSV24BE		PSV24DCBE	

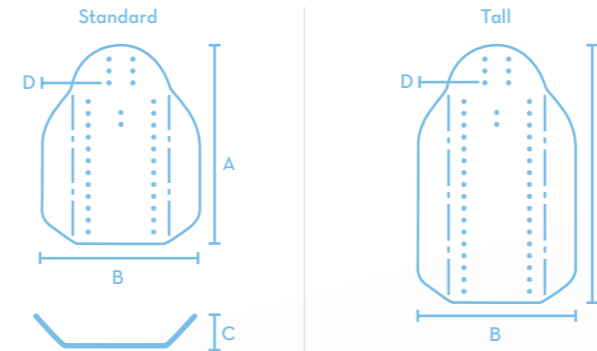

Hole spacing is 1 inch on centre starting from the bottom up

STANDARD SHELL DIMENSIONS

SHELL SIZE	A - HEIGHT	B - WIDTH	C - DEPTH	D - NO. OF PRE-DRILLED HOLES
30.5cm (12")	31.1cm (12.25")	23.5cm (9.25")	5.1cm (2")	9
35.6cm (14")	37.5cm (14.75")	28.6cm (11.25")	5.7cm (2.25")	11
40.6cm (16")	42.5cm (16.75")	33.7cm (13.25")	6.4cm (2.5")	12
45.7cm (18")	46.4cm (18.25")	38.7cm (15.25")	7cm (2.75")	13
50.8cm (20")	47.6cm (18.75")	43.8cm (17.25")	8.3cm (3.25")	13
55.8cm (22")	50.8cm (20")	48.9cm (19.25")	8.9cm (3.5")	13
61cm (24")	53.3cm (21")	54cm (21.25")	9.2cm (3.625")	14


TALL SHELL DIMENSIONS

SHELL SIZE	A - HEIGHT	B - WIDTH	C - DEPTH	D - NO. OF PRE-DRILLED HOLES
30.5cm (12")	38.7cm (15.25")	23.5cm (9.25")	5.1cm (2")	12
35.6cm (14")	47.6cm (18.75")	28.6cm (11.25")	5.7cm (2.25")	15
40.6cm (16")	55.2cm (21.75")	33.7cm (13.25")	6.4cm (2.5")	17
45.7cm (18")	59.1cm (23.25")	38.7cm (15.25")	7cm (2.75")	18
50.8cm (20")	59.1cm (23.25")	43.8cm (17.25")	8.3cm (3.25")	18



Leaflet Code: LL223

For more information contact your local dealer:



driveTM
www.drivemedical.co.uk

All sizes and weights stated are nominal. All errors and omissions excepted.
Due to the limitations of photography and the printing process colours shown may not be 100% accurate.