

Invacare® | **Flo-tech® Contour Visco**



3 Year warranty



Conforms to the shape of the body

The *Invacare Flo-tech Contour Visco* cushion is moulded from visco-elastic (memory) foam, which is temperature-sensitive to the body. As the foam warms, the contoured surface area conforms to body shape, improving posture, increasing comfort and maximising the surface area in contact with the client, thereby optimising pressure reduction. The *Flo-tech Contour Visco* is suitable for those considered at 'High Risk'* of developing pressure ulcers.



Yes, you can.®

Deep moulded visco elastic (memory) foam cushion

The temperature-sensitive foam conforms to the client's body, maximising the contact area, delivering exceptional levels of comfort and optimising pressure reduction. The deep foam makes this cushion particularly suitable for heavier and taller clients.



Contoured seating surface withommel and side support

Offers increased lateral support and pelvic stability to enhance posture.

Features and Options



Water Resistant Cover

Cover is multi stretch, water-resistant and vapour-permeable. Cover features an anti-slip base and carry handle.

Available in a wide range of sizes

To view the cushion size range, please visit www.invacare.co.uk.

Four foam densities available

Available in four densities to accommodate a range of client weights -
 Soft: 8 stone (up to 51 kg),
 Medium: 12 stone (up to 76 kg),
 Firm: 16 stone (up to 102 kg) and
 Very Firm: 20 stone (up to 125 kg).

Technical data



Flo-tech Contour Visco A range of 51 sizes



Height: 10.3cm / 4"



Weight: 1.60 kg



Max. User Weight: 125 kg (20 stone)



Cleaning: 80°


Water Resistant
 Vapour Permeable

Manufactured to European Standards:
 Interior Foam: BS EN 597-1, BS EN 597-2, BS 5852: 1982 Crib 5
 Cover: BS 7175: 1989 Crib 5


Colours



Black Smooth

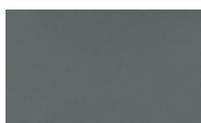

Black Wipedown

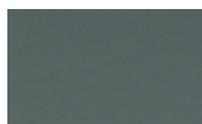

Black Towelling



Black Spacer Fabric



Kingfisher Green

Warranty 3 Years

INVACARE® Limited
 Pencoed Technology Park - Pencoed -
 Bridgend CF35 5HZ
 Tel: +44 1 656 776 222
 Fax: +44 1656 776 220
 E-mail: uk@invacare.com
www.invacare.co.uk
 Sales Order E-mail:
ordersuk@invacare.com

INVACARE® Ireland Ltd
 Unit 5 - Seatown Business Campus
 Seatown Road - SWORDS - County
 Dublin
 Tel. +353 1 8107084
 Fax +353 1 8107085
 Email: Ireland@invacare.com
www.invacare.ie

© 2010 Invacare® International Sarl
 All rights reserved.

All information quoted is believed to be correct at time of print. Invacare® reserves the right to alter product specifications without prior consultation.

Flo-tech Contour Visco
 UK - 10/2010 - 1505613

* Essential nursing care is pivotal in pressure ulcer prevention. This product will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Yes, you can.®