

Invacare® | ***Flo-tech® Contour***



3 Year warranty




Designed for comfort and improved posture

The surface of the *Invacare Flo-tech Contour* cushion is contoured into sections. This unique design improves posture and increases the support surface area, effectively dispersing weight away from specific pressure points. The combination of improved posture and increased support surface area will contribute significantly to pressure ulcer prevention. This deep cushion is suitable for heavier clients (up to 30 stone), and those considered 'At Risk'* of developing pressure ulcers.



Yes, you can.®

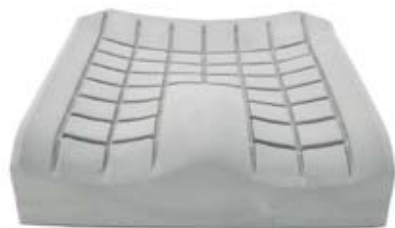

Moulded foam cushion

Moulding allows **Flo-tech** cushions to provide a 3D form, which cannot be emulated in cut foam. Moulded foam is highly durable, increasing the product's longevity. The pommel at the front of the cushion separates the thighs, encouraging good posture and allows air to circulate.

Contoured seating surface

Offers improved stability and postural support. The castellations of the cushion allow slight movement which in turn reduces pressure. The contoured surface helps to prevent the client sliding forward, therefore reducing the potential for shear and friction forces.

Features and Options



Water Resistant Cover

Cover is multi stretch, water-resistant and vapour-permeable. Cover features an anti-slip base.

Available in a wide range of sizes

To view the cushion size range, please visit www.invacare.co.uk.

Technical data

Flo-tech Contour



A range of 83 sizes



Height:
10.3 cm / 4"



Weight:
1.38 kg



Max. User Weight:
191kg (30 stone)



Cleaning: 80°

Manufactured to European Standards:
BS EN 597-1 Cig, BS EN 597-2 Match, and BS 7176
(Crib 5)

Warranty 3 Years

Colours



Black Smooth



Black Wipedown



Black Towelling



Black Spacer Fabric



Kingfisher Green

INVACARE® Limited
Pencoed Technology Park - Pencoed -
Bridgend CF35 5HZ
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
www.invacare.co.uk
Sales Order E-mail:
ordersuk@invacare.com

INVACARE® Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County
Dublin
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2010 Invacare® International Sarl
All rights reserved.

All information quoted is believed to be correct at time of print.
Invacare® reserves the right to alter product specifications without prior consultation.

Flo-tech Contour - UK - 10/2010
1505440



Yes, you can.®

* Essential nursing care is pivotal in pressure ulcer prevention. This product will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.