


Invacare® | **Flo-tech® Solution**



**3 Year warranty**

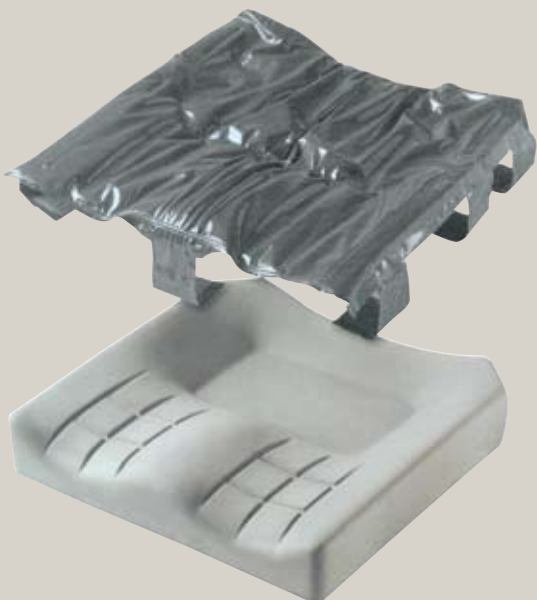

**Superior protection for 'Very High Risk' clients**


The *Invacare Flo-tech Solution* cushion uses a combination of a slim, contoured foam base with a dual layer of gel sacs over the entire cushion surface. The dual gel sacs and the unique moulded foam shape of the *Flo-tech Solution* significantly increases the support surface area. This delivers a superior level of protection for those considered to be at 'Very High Risk'\* of developing pressure ulcers.



**Yes, you can.®**

**Delivers superior protection from pressure ulcers**

The *Invacare Flo-tech Solution's* design increases the contact area between the client and the foam cushion. The gel re-distributes pressure away from bony prominences, offering increased protection. The additional layer of gel is incorporated under the ischial and sacral area to further enhance the pressure reducing qualities of this cushion, making this cushion suitable for 'Very High Risk' clients.

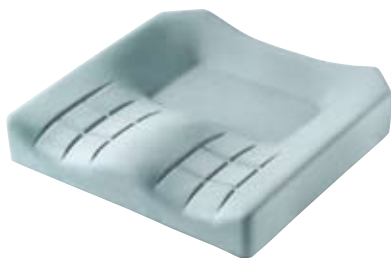

**Slimline contoured cushion**

Some clients prefer a slimline cushion, finding it more compatible with their daily living activities. The low profile is particularly beneficial for the facilitation of side transfers. This may make this cushion more suitable for use with existing equipment.

**Features and Options**



**Water Resistant Cover**

Cover is multi stretch, water-resistant and vapour-permeable. Cover features an anti-slip base and carry handle.

**Available in a wide range of sizes**

To view the cushion size range, please visit [www.invacare.co.uk](http://www.invacare.co.uk).

**Invacare Flo-tech Deep Solution (Option)**

Option of a 4" (10 cm) "Deep" foam base is also available in a range of sizes.

**Technical data**



A range of 44 sizes



Height: 7.8 cm and 10.2cm (3" Standard, 4" Option)



Weight: 3.08 kg



Max. User Weight: 140 kg (22 stone)



Cleaning: 80°

**Flo-tech Solution**

Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match, and BS 7176 (Crib 5)

INVACARE® Limited  
Pencoed Technology Park - Pencoed - Bridgend CF35 5HZ  
Tel: +44 1 656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)  
Sales Order E-mail: [ordersuk@invacare.com](mailto:ordersuk@invacare.com)

**Colours**

**Warranty 3 Years**



Black Smooth



Black Wipedown



Black Towelling



Black Spacer Fabric



Terracotta Swatch (Only available for welded wipe down covers)



Teal blue PU

INVACARE® Ireland Ltd  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS - County Dublin  
Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [Ireland@invacare.com](mailto:Ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

© 2010 Invacare® International Sarl  
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare® reserves the right to alter product specifications without prior consultation.

Flo-tech Solution - UK - 10/2010  
1505703

\* Essential nursing care is pivotal in pressure ulcer prevention. This product will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



**Yes, you can.®**