



Stability and safety meet performance in style

The remarkable *Invacare TDX SP2^{NB}*, with its compact narrow base, boasts durability and ultimate performance. Great manoeuvrability indoors and out, with a wide choice of seating for optimum comfort and support makes the *TDX SP2^{NB}* the powerchair for individuals who need full functionality, but do not want to compromise on style. Key features include:

- Enhanced stability
- Intuitive driving experience
- Simply Smart LiNX control system
- Total Drive Experience



LiNX®



Yes, you can.®



Compact and Versatile

The **TDX SP2^{NB}** is designed for ultimate freedom. Whether effortlessly manoeuvring through tight spaces or tackling challenging outdoor terrain, performance and comfort are guaranteed.

The Perfect Combination

Blending the most high performance centre wheel drive base on the market with the most highly adjustable and adaptable seating system, the **TDX SP2^{NB}** is completed with the new **LiNX** control system. Stylish, reliable and intuitive.

Simply SMART LiNX Control System

Invacare **LiNX** is our new insight inspired control system that learns over time to bring a revolutionary driving experience. Thanks to its unique modular design, this futureproof system can be seamlessly expanded and tailored to complement changing needs.

This simply smart technology allows programming and maintenance updates to be made wirelessly via bluetooth and in real time, making it easy to service. With easy access to statistical and diagnostic information, **LiNX** will allow complete visibility. Features include:



Enhanced interface

The design of the controls has focused on ease of use, simplicity and reliability. The Simply SMART controls are easy to understand and use for everyone, including those with limited cognitive abilities. All controls feature an accurate battery gauge, a lower force joystick, a physical speed pot and an icon based menu for controls with lights and powered positioning functionality.



Direct access to seating functions

This is easily accessible via the optional 10-way switch and suitable for those who have difficulty reaching the buttons on the control or with their hand positioning.



Total control

For more complex control needs, **LiNX** improves functionality, driving control and performance while **G-Trac[®]** (optional) ensures the directions indicated by the control are followed by the chair.

FEATURES & OPTIONS



SureStep® Suspension

SureStep suspension means all six wheels of the **TDX SP2^{NB}** stay firmly on the ground and the seat stays level, no matter the terrain.



Stability Lock

Stability Lock comes into play when base and front castors start tipping forward in relation to the rear castors, providing stability and safety.



Dahl Docking Station

The **TDX SP2^{NB}** can be prepared to fit the small and easy to use Dahl docking station enabling even simpler and safer transportation.



Powered fixed pivot tilt option

The simplified powered tilt option allows the seat angle to be repositioned by up to 20°.



Flip-up armrest

The **Modulite™** seat can be combined with a multi-adjustable flip-up armrest for personalised comfort and easier transfers.



Side rails

Side rails on each side of the chair are mounted to telescopic bars and allow for the individual placement of armrests and accessories.



Seat plate settings for width and depth

The 3-piece seat plate covers all possible width and depth settings.



Compact Turning Radius

With a compact turning radius of 560mm, the **TDX SP2^{NB}** is great indoors or out.



Telescopic seat frame

The telescopic seat frame is at the heart of the **Modulite™** seat design and enables a whole range of positioning options.



Powered tilt-in-space

A maximum seat angle of 30° ensures comfort and safety as well as including a centre of gravity shift.

1/ Invacare **Modulite™** seat components are cross compatible across the Invacare power wheelchair range and partly with the Invacare Rea® family.

Technical data

Invacare TDX SP2^{NB}	380 - 530 mm ²	410 - 510 mm 510 - 560 mm ³	450 / 475 mm ^{4,5} 475 - 775 mm ⁶	480, 540, 560 mm ⁷	245 - 310 mm	290 - 460 mm ⁸	0 - 20° / 0 - 30° ⁵
Invacare TDX SP2^{NB}	90° - 120°	1010 - 1170 mm ¹¹	610 mm ¹⁰	1160 mm ⁹	980 mm	150 kg	50 Ah
Invacare TDX SP2^{NB}	1120 mm	75 mm	9° / 16 %	2 x 340 W	up to 26 km ¹²	6 and 10 km/h	138 kg ¹¹


- 1/ Seat cushion 430 mm
- 2/ Additional 50 mm by armrest adjustment
- 3/ Extended profile rail package for larger seat depth
- 4/ Ground to seat plate without cushion
- 5/ With seat lifter or powered tilt
- 6/ With lifter
- 7/ Measured from seat plate
- 8/ Measured from footplate to seat plate
- 9/ Depending on seat position
- 10/ Width of wheelbase
- 11/ Depending on configuration

- 12/ According to ISO 7176-4:2008
- 13/ Invacare **Modulite[™]** seat components are cross compatible across the Invacare power wheelchair range and partly with the Invacare Rea[®] family.
- 14/ Compliant with the directive on the Restriction of the use of certain Hazardous Substances in electrical and electronic equipment - RoHS (2011/65/EC).
- 15/ Compliant with the regulation concerning the Registration, Evaluation, Authorization and Restriction of Chemicals.



matrix[®] SEATING SERIES Invacare **Modulite[™]**

Shroud colours

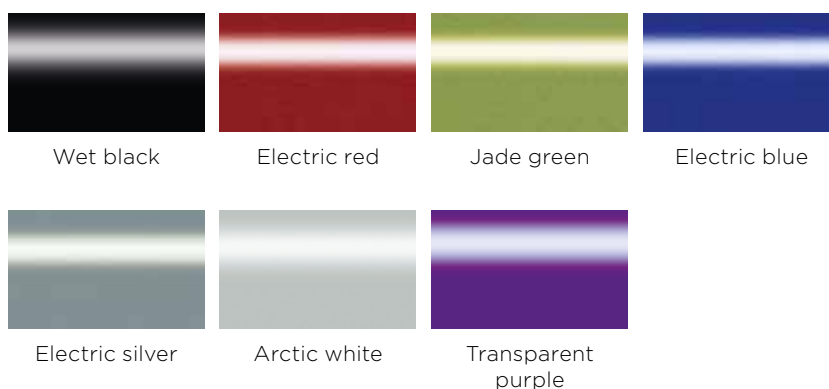

The LiNX word, mark and logos are registered trademarks owned by Dynamic Controls and any use of such marks by Invacare is under license.

For more information about this wheelchair, please see www.invacare.co.uk. For information about typical description, intended use, adjustments and transportation, please refer to the user guide.

Invacare Limited
Pencoed Technology Park -
Pencoed - Bridgend CF35 5AQ -
United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS -
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2016 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. The images in this brochure are for illustration purposes only and may vary slightly from the exact product available.

TDX SP NB2 - UK - 06/2016 -
1578272



Follow us on Twitter:
[@Think_Mobility](https://twitter.com/Think_Mobility)



Follow us on YouTube:
Invacare Ltd - UK



Like us on Facebook:
Invacare UK



Yes, you can.[®]