



### A versatile compact chair for people with changing needs

The *Invacare*® Rea **Dahlia** has it all, a slim and compact wheelbase in a modern design. It incorporates an adjustable centre of gravity option and a light and easy tilt and recline function. For those seeking mobility with comfort and pressure relief, this versatile chair is perfect for both passive and active clients alike.

The **Dahlia** features two tilt angles of either 30° or 45°, which improves the positioning capabilities as the tilt and recline options reduce pressure. Component Cross-Compatibility (CCC) across the Rea product range includes the **Flex3** seating.





### Weight shifting with ease

The 45° option provides greater tilt possibilities, delivering excellent pressure relieving properties. The unique construction of the tilt mechanism makes it light to adjust, requiring little muscle force.


### An ideal chair for foot propelling

The 30° tilt option is suitable for more able clients who can readily propel themselves. The seat height can be set as low as 325mm, allowing the possibility to foot propel. There is also a self tilt option which allows an individual the ability to change their posture to suit themselves.

### Benefits for both *Dahlia* 30° and 45°:

- Easy to manoeuvre - self propel rear wheels, easily accessible
- Small overall width - narrow footprint
- Centre of gravity adjustment allows better weight distribution, allowing lighter pushing and easy turning capabilities.
- Stable - DSS, Dual Stability system maintains chair stability in every tilt position
- Light and smooth tilt and recline action
- Component Cross-Compatibility from the Invacare® family

### Features and Options

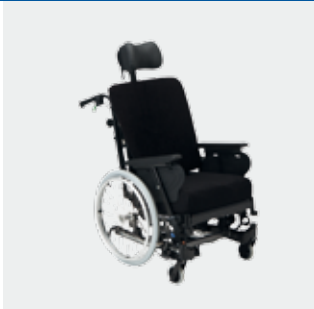

**Dahlia** 30° can be foot propelled with seat heights as low as 325mm.



**Dahlia** 30° with self tilt for the more able person to actively change posture.



**Dahlia** 45° for optimum pressure relief and weight distribution.



The Dual Stability System, DSS helps to avoid friction and shear when reclining.



## Features and Options




On the 30° version of the **Dahlia**, in the first 10° of tilt the front seat height is only raised a small amount, this allows the feet to remain close to the floor.



Moulded seats, or individually tailored seating is possible.



The **Matrx** is available as an option. It provides a stable and firm support to the body. The **Matrx Elite Deep** has a built-in lateral support and the **Matrx Elite** has a mild contour, which allows a greater range of movement.

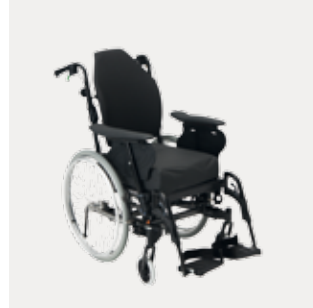

For a wider backrest on a smaller seat frame, backrest joint width extensions are an option offering an extra 50mm of width on the seat.



Seat plate extensions offer individual support where needed.



Adjustable seating independent of depth and width.



**Dahlia** has an optional tension adjustable backrest.



The overall width of the **Dahlia** is comparable to an active chair. This makes it easy to manoeuvre in confined spaces and through narrow doorways.



Kerb lever is always mounted on the right of the chair.



The adjustable centre of gravity on the rear wheels allows weight distribution for easier and lighter pushing.



Anti tippers when not required can be flipped up. When in use they rest at a 5° outward angle. This allows the attendant to walk closer to the chair without colliding with the anti tipper.

## Technical data



Rea **Dahlia** 390-450/440-500/  
490-550 mm  
16" (390-490/440-540/  
490-590 mm)



430 - 500 mm



30° : 330 - 500 mm  
45° : 350 - 500 mm



500 - 800 mm  
(without seat  
cushion)



250 - 320 mm



330 - 500 mm



30°: -1 - 30°  
45°: 0 - 45°



Rea **Dahlia**

0° - 30°



SW + 160 mm (16")  
SW + 180 mm (24")



30°:  
825 - 1250 mm  
45°:  
850 - 1250 mm



900 - 1180 mm



>35kg<sup>1</sup>



Max 135 kg



Transport weight  
21/25 kg (SW 39 cm)



Rea **Dahlia**

755 mm




Crash-tested  
ISO 7176-19

The Rea **Dahlia** complies with ISO 7176-19 specifications. Note that special conditions apply regarding accessories. For detailed information regarding vehicular transportation and for more pre-sale information for this wheelchair, please consult the following pages on our website: [www.invacare-rea.com](http://www.invacare-rea.com) under Rea wheelchairs/Technical downloads.

For options and accessories, see "Rea Accessory" brochure.

<sup>1/</sup> Overall weight in seat width 390mm and with lightest configuration

## Frame colour

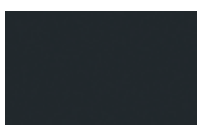

Black

## Upholstery



Black Dartex  
TR26



Black Spacer  
TR36

<sup>2/</sup> Component Cross-Compatibility with the Invacare® Rea family

Rea, Rea design and DSS (Dual Stability System) design are registered trademarks of Invacare International.

The Rea **Dahlia** is designed with DSS (Dual Stability system). It provides the user with enhanced stability and safety when being tilted or the backrest is reclined.



Invacare Limited  
Pencoed Technology Park -  
Pencoed - Bridgend CF35 5AQ -  
United Kingdom  
Tel: +44 1656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail:  
[ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Ireland Ltd  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS -  
County Dublin - Ireland  
Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [ireland@invacare.com](mailto:ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

© 2016 Invacare International Sàrl  
All rights reserved.

All information quoted is believed to be correct at time of print.  
Invacare reserves the right to alter product specifications without prior consultation.

Dahlia - UK - 02/2016 - 1576793