

Invacare® | **Propad® Premier Cushion**



3 Year warranty



Designed for comfort, support and pressure reduction

The *Invacare Propad Premier Cushion* has been designed to give comfort and pressure reduction. The cushion has a dual layer construction to ensure improved support as well as prolonged longevity. The **Propad Premier Cushion** is suitable for those considered 'At Risk'* of developing pressure ulcers.




Yes, you can.®



Features and Options

Dual-layer profiled foam cushion, designed to deliver enhanced pressure care

The lower layer of the *Invacare Propad Premier Cushion* is constructed from a resilient high density foam and is shaped to deliver increased stability and support for the client. The upper section of the cushion features a layer of castellated foam. Each castellation moves independently to deliver enhanced support and pressure area care to the client, as well as ensuring good air circulation, reducing the potential for sweating and associated skin maceration.

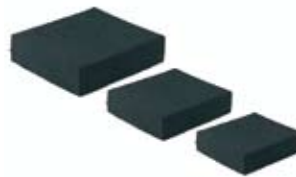

Cover helps to prevent a build up of shear and friction forces

Cover is multi stretch, water-resistant and vapour-permeable. Cover features an anti-slip base.

Available in a wide range of sizes

To view the cushion size range, please visit www.invacare.co.uk.

Technical data



Wide range of sizes



Height: 102 mm / 4"



Weight: 1 kg



Max User Weight: 110kg / 17.3 stone



Cleaning: 80°

Propad Premier Cushion

Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match, and BS 7175 (Crib 5)

Warranty 3 Years

INVACARE® Limited
Pencoed Technology Park - Pencoed - Bridgend CF35 5HZ
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
www.invacare.co.uk
Sales Order E-mail: ordersuk@invacare.com

INVACARE® Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County Dublin
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2010 Invacare International Sarl
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

Propad Premier Cushion
UK - 10/2010 - 1505221



Yes, you can.®

Colours



Blue PU



Black Towelling



Black Smooth

* Essential nursing care is pivotal in pressure ulcer prevention. This product will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

