

rea | Azalea Family



Body length compensation



Weight shifting



One big family

Invacare Rea Passive chairs offer people the perfect fit and comfort, with a wide selection of models, options and compatibility across the range.

Allowing for customers who have specific needs, *Rea* has developed a range of niche models, each chair highly adjustable to provide a better fit for the individual.

All chairs can be quickly adjusted, and come with the *Invacare Dual Stability System*® (DSS) and are proven to be safe and durable to last.













Yes, you can.®



In order to provide the client with the best fit, the **Azalea** offers individual size adjustments to improve posture. In addition to the popular width and depth adjustment of the seat area, we also offer the same adjustment on the backrest. This provides a unique combination for postural support and comfort as the adjustments are independent of one another.

Define the posture support you are looking for and it can be fulfilled on the **Rea Azalea Family**.

Suggestion of backrest options, backrest plate or tension adjustable

Plate (Tension adj)	Plate (Tension adj)	Plate Tension adj	Plate (Tension adj)	(Plate) Tension adj	(Plate) Tension adj	(Plate) Tension adj	Tension adj	Tension adj	Tension adj	(Plate) Tension adj
										
Short	Tall	Average	Wide	Wide above	Wide below	Larger in front	Larger behind	Kyphotic back	S-shaped back	
All fittings as close as possible	Extended backrest height and seat depth	Fittings of backrest and seat in standard positions	Backrest and seat as wide as possible	Upper part of backrest is wide apart, lower seat is narrow	Lower part of backrest is wider as the seat	Upper part of backrest turned to allocate space around shoulders	Lower part of backrest turned to allocate space for the protruding bottom	Midsection of the backrest is opened up with a tighter support in the lower regions	Low midsection is tight whereas the other parts follows the S-shaped spine	

Backrest options



The Flex3 tension adjustable option can be adjusted in both width and height. In addition the backrest posts can be angled and the whole backrest can be turned upside down to create space where it's needed most.



The Flex3 backrest plate option is height and width adjustable.



Features and options

Offers excellent seating comfort and many adjustment possibilities!

The **Rea Azalea** is the natural choice for the moderately active client as well as for the more passive. The comfortable seating system allows you to alter the seat position, ensuring maximum comfort at all times. This combined with good driving characteristics makes **Azalea** the ideal choice.

Tilt & recline makes life comfortable

The **Rea Azalea** has a 25° seat tilt and a 30° backrest recline. Both tilt and recline can be easily adjusted manually or electrically by the assistant and also by the client on the electrical version.



Rea Azalea Self Propel

The original self propelling wheelchair with large rear wheels.



Rea Azalea Transit

The transit shares the frame with the **Rea Azalea** Self Propel but is mounted with small rear wheels.

Technical data

Rea Azalea SP Rea Azalea TR	390-450/440-500/ 490-550 mm (390-490/440-540/ 490-590 mm)	430 - 500 mm	400/450 mm (450 mm)	Backrest plate: 560-790 mm Tension adjustable backrest: 540-715 mm	240 - 360 mm	330 - 500 mm	-1° to 23° (+1° to 25°)
Rea Azalea SP ReaAzalea TR	0°-30°	SW + 220 mm (SW + 230 mm)	900-1250 mm	950-1020 mm	34 kg (SW 440 mm)	Max 135 kg Transport weight: 16,5 kg	
	The Rea Azalea and Azalea Assist complies with ISO 7176-19 specifications. Note that special conditions apply regarding accessories. For detailed information regarding vehicular transportation and for more pre-sale information for this wheelchair, please consult the following pages on our website: www.invacare-rea.com under Rea wheelchairs/Technical downloads.						
Rea Azalea SP ReaAzalea TR	Crash-tested ISO 7176-19						

Frame colour



Grey



Azurite Blue

Upholstery



Black Dartex
TR26



Black Spacer
TR36



Features and options

Rea Azalea Base

Designed for people seeking a custom made seat combined with a first-class wheelchair base, the **Rea Azalea Base** is adapted from our standard **Rea Azalea** model and boasts all the advantages of this sturdy, tilting chair, such as the special weight-shifting mechanism. A key feature of the **Base** is the 'Tilt in Space' function - a weight-shifting mechanism that ensures stability for the client when tilted.

Rea Azalea and Alber viaplus V12

Unlike many other pushing aids, the **viaplus V12** does not impede the adjustability of a tilt in space wheelchair (with a seat height of 450 mm). All of the wheelchair's functions remain accessible. That's why the **viaplus V12** can simply be left on the wheelchair without restricting the use of its functions.



Rea Azalea Base

Ideal for those in need of a stable, secure chair, with custom-made comfort. The **Base** can be used together with all parts and accessories of other **Rea Azalea** models (except trunk support).



Rea Azalea & Alber viaplus V12

Technical data

Rea Azalea Base	340*/390/440/ 490 mm (+ 50/100 mm)	380 - 450 mm* 430 - 500 mm	400/450 mm	600* / 700 mm	240 - 360 mm	330 - 500 mm
Alber viaplus V12				600* / 700 mm	240 - 360 mm	330 - 500 mm
						-1° to +23° (+1° to 25°)
Rea Azalea Base	0°-30°	SW + 220 mm (SW + 230 mm)	1300* / 1400 mm	930 - 1040* / 950 - 1070 mm	20 kg	Max 135 kg
Alber viaplus V12	0°-30°	SW + 220 mm (SW + 230 mm)	1300* / 1400 mm	900 - 1020* / 950 - 1070 mm	28 kg* (SW 340) / 34 kg (SW 440 mm)	Max 75*
						Transport weight: 15 kg
						15* / 18,7 kg



Crash-tested
ISO 7176-19

* **Rea Azalea Minor** has smaller dimensions.
The **Rea Azalea Base** complies with ISO 7176-19 specifications. Note that special conditions apply regarding accessories. For detailed information regarding vehicular transportation and for more pre-sale information for this wheelchair, please consult the following pages on our website: www.invacare-rea.com under Rea wheelchairs/ Technical downloads.

Frame colour



Grey



Azurite Blue

Upholstery



Black Dartex
TR26



Black Spacer
TR36



Features and options

Rea Azalea Minor – Designed for petite adults and teenagers

The **Rea Azalea Minor** is designed for people who require a seat that has a smaller depth and width. Adapted from the **Rea Azalea**, the **Azalea Minor** boasts all the advantages of a reliable, tilt-able wheelchair with its special weight-shifting mechanism. It is also available with a tension adjustable backrest.

Rea Azalea Tall - No limits on comfort for taller person

Rea Azalea Tall is specially designed to meet the needs of tall clients who require a Tilt in Space wheelchair with a longer seat support. Adapted from the **Rea Azalea**, the **Rea Azalea Tall** boasts all the advantages of a reliable, tilting wheelchair and offers a unique weight shifting mechanism.



Rea Azalea Minor
The push handles of the **Azalea Minor** are angled outwards, increasing comfort for carers.



Rea Azalea Tall
The **Rea Azalea Tall** at the top of its class for comfort and stability.

Technical data

Rea Azalea Tall	390-450/440-500/ 490-550 mm (390-490/440-540/ 490-590 mm)	480 - 550 mm	500 mm	Backrest plate: 560 - 790 mm Tension adjustable backrest: 540 - 715 mm	240 - 360 mm	330 - 500 mm	-1° to 23° (+1° to 25°)
Rea Azalea Minor	340*/390/440 / 490 mm (+50/100 mm)	380 - 450 mm* 430 - 500 mm	400/450 mm	550 - 650 mm	240 - 360 mm	330 - 500 mm	+1° to 25°
Rea Azalea Tall	0°-30°	SW + 220 mm (SW + 230 mm)	1000 - 1300 mm	950 - 1070 mm	34 kg (SW 440 mm)	Max 135 kg	Transport weight 18,7 kg
Rea Azalea Minor	0°-30°	SW + 220 mm (SW + 230 mm)	900 - 1250 mm	900-970mm	32 kg	Max 75* / 135 kg	14,5 kg



Crash-tested
ISO 7176-19

* **Rea Azalea Minor**
Rea Azalea Tall and Azalea Minor complies with ISO 7176-19 specifications. Note that special conditions apply regarding accessories. For detailed information regarding vehicular transportation and for more pre-sale information for this wheelchair, please consult the following pages on our website: www.invacare-rea.com under Rea wheelchairs/Technical downloads.

Frame colour



Grey



Azurite Blue

Upholstery



Black Dartex
TR26



Black Spacer
TR36



Features and options

A tailor-made solution for heavier clients

The **Rea Azalea Max** is designed to meet the specific needs of the heavier person. The **Azalea Max** is a niche product and a complementary addition to the existing **Azalea family** from Invacare.

Continued Compatibility

The **Rea Azalea Max** offers the same features as the rest of the **Rea Azalea family** such as tilt-in-space with weight-shifting and DSS[®] (Dual Stability System) mechanism. Stability in this wheelchair originates from body length compensation in the back and leg rests, beneficial for people with limited movement and have difficulty changing position.



Designed to meet the specific needs of people with extra weight.



Two-way backrest

The lower and upper parts of the shaped backrest are interchangeable, giving the option of more space in the seat or chest area.



Power tilt

The hand control for tilt and recline can be used by the attendant or the person using the chair.

Technical data

Azalea Max	550; 610 + 50 mm	530 - 600 mm	400 / 450 mm	500 + 80 mm	320 - 420 mm	330 - 500 mm	1° - 20°
Azalea Max Assist	550; 610 + 100 mm	530 - 600 mm	400 / 450 mm	500 + 80 mm	320 - 420 mm	330 - 500 mm	1° - 20°
Azalea Max	0°-30°	SW + 250 mm	1020 - 1240 mm	1050 - 1120 mm	54 kg		
Azalea Max Assist	0°-30°	SW + 220 mm	1020 - 1240 mm	1050 - 1120 mm	54 kg		
		<p>For other Rea Azalea accessories, please consult the Rea Azalea brochure. The Rea Azalea Max complies with ISO 7176-19 specifications. Note that special conditions apply regarding accessories. For detailed information regarding vehicular transportation and for more pre-sale information for this wheelchair, please consult the following pages on our website: www.invacare-rea.com under Rea wheelchairs/ Technical downloads.</p>					
Azalea Max	33 kg	max. 180 kg					
Azalea Max Assist	37 kg	max. 180 kg					

Invacare Limited
 Pencoed Technology Park - Pencoed - Bridgend CF35 5AQ - United Kingdom
 Tel: +44 1 656 776 222
 Fax: +44 1656 776 220
 E-mail: uk@invacare.com
 Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk


Invacare Ireland Ltd
 Unit 5 - Seatown Business Campus
 Seatown Road - SWORDS - County Dublin - Ireland
 Tel. +353 1 8107084
 Fax +353 1 8107085
 Email: Ireland@invacare.com
www.invacare.ie

© 2015 Invacare International Sarl
 All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

Rea Azalea Family - UK - 07/2015 - 1554736

Frame colour

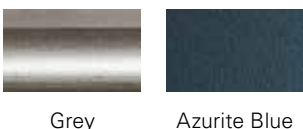

Grey Azurite Blue

Upholstery



Black Dartex TR26 Plush

Rea, Rea design and DSS (Dual Stability System) design are registered trademarks of Invacare International. The **Rea Azalea** is designed with DSS (Dual Stability system). It provides the user with enhanced stability and safety when being tilted or the backrest is reclined.



Yes, you can.[®]



Follow us on Twitter:
[@Think_Mobility](https://twitter.com/Think_Mobility)



Follow us on YouTube:
 Invacare Ltd - Think Mobility